

Mobi-chair™ Floating Beach Wheelchair











Although the Mobi-chair™ floats, it is not a boat. The sitting position is stable and you won't slip out of the Mobi-chair ™ under normal conditions of use.

However, for safety reasons, a life jacket is required for non swimmers (this may also be required by law). Mobi-chair ™ is not intended for areas with current or strong waves.



ADVANTAGES



UV and sea water resistant

- > Mobile
- > Safe & robust
- Floating
- > Anticorrosion
- > Useful & ergonomic

- > For both adults and children
- > Seamless transition from boardwalk-to beach-to water

Technical specificities

| Folding frame | Aluminum frame & stainless fork and hardware, salt water resistant | | | |
|------------------------|--|---------------|---------------|-------------------|
| Mobi-chair dimensions® | Width | Length | Height | Weight |
| | 107 cm 42" | 152 cm 60" | 127 cm 50" | 30kg 66 lbs |
| Seat dimensions | Width | Depth | Height | Max load capacity |
| | 46 cm 18" | 43 cm 17" | 48 cm 19" | 136 kg 300 lbs |
| Adjustable backrest | 3 positions: relaxation, seated, swimming | | | |
| Seat color | Blue | | | |
| Seat Fabric | UV Resistant - Water Resistant Ventilated - Non Allergenic | | | |
| Detachable Armrest | PVC Flip up - detachable floating armrest - low pressure | | | |
| Front wheel | 30cm x 16cm 12" x 6-1/2» Rubber Pneumatic Tire*1 w/quick release | | | |
| Rear wheel | 40cm x 19cm 16" x 7-1/2» Rubber Pneumatic Tire*2 w/quick release | | | |

Adaptable & high quality wheelchair




- Easy mobility through sand and other rough surfaces
- Safety harness



Rubber pneumatic tires absorb the shock and discomfort caused by rugged terrain



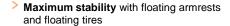

> Handles integrated into the armrests



Transportation and storage are

provide support

convenient

Quick installation & low maintenance

QUICK INSTALLATION


Assembly and disassembly of Mobi-chair™ without material and in 10 minutes maximum.

Completely foldable, Mobi-chair™ is easily transportable and fits perfectly in the trunk of a car.

EASY **M**AINTENANCE

Once a week, lubricate all connecting rods and safety screws.

After each use, rinse Mobi-chair® with fresh water to avoid corrosion and store it in a clean, dry place.

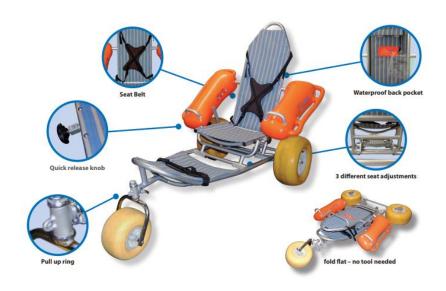







Multiple applications

SAND

OCEANS & SEAS





LAKES



SWIMMING **P**OOLS









